

Distributed By:

**MEDI-STIM** INC.

**VIVATENS**  
OTC

Prescription T.E.N.S.  
NOW AVAILABLE  
**OVER-THE-COUNTER**  
FOR DRUG-FREE PAIN RELIEF!

**VIVATENS**  
*Mini* OTC




Digital microprocessor controlled devices with a lumbar belt that has been cleared by FDA to provide temporary relief of pain associated with sore or aching muscles of the lower back due to strain from exercise or normal household and work activities.

No larger than a cell phone, these portable battery-operated devices deliver intermittent pulses that can help block pain signals and stimulate the body's production of endorphins, a natural pain killer released by the brain.

Stimulation from the T.E.N.S. device is delivered through a lumbar belt with a built-in 2 channel cable system and electrode pads. 8 preset programs designed for professional pain relief offers effective and easy treatment.

When used properly, the VivaTENS OTC Back Pain Relief System can provide pain relief lasting for several hours. TENS treatment can be combined with other therapy, such as exercise, to achieve optimum health benefits.

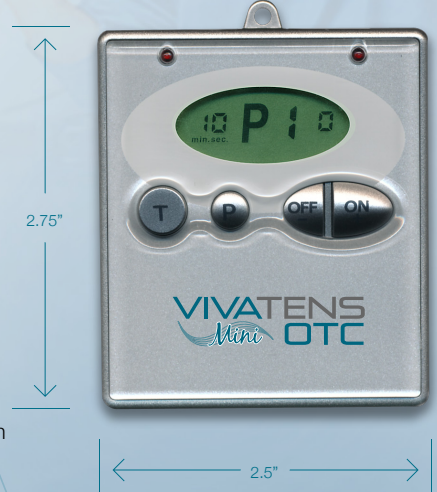

- ▶ Dual Channel
- ▶ 8 preset programs
- ▶ Fully Adjustable Intensities
- ▶ Treatment Timer: 10, 20, 30, 60 min. or Continuous
- ▶ LCD Channel Power Indicators
- ▶ Power Supply: 9V Battery
- ▶ Size: 3.75"(L) x 2.55"(W) x 1.4"(D)
- ▶ Weight: 150 grams with battery
- ▶ 1-year warranty



- ▶ Dual Channel
- ▶ 8 Preset Programs
- ▶ 10 Intensity Levels
- ▶ Treatment Timer: 10, 20, 30, 60 Minutes
- ▶ LED Channel Power Indicators
- ▶ Power Supply: 2 ea AAA Alkaline Batteries
- ▶ Size: 2.5"(L) X 2.75"(W) X 0.6"(D)
- ▶ Weight: 60 Grams (Including Batteries)
- ▶ 1-Year Warranty



	Intensity	Pulse Width	Pulse Rate	Mode
P1	0~50 mA	260 uS	15 Hz	Constant
P2	0~50 mA	260 uS	100 Hz	Burst
P3	0~50 mA	260 uS	120 Hz	Constant
P4	0~50 mA	150-260 uS	2-100 Hz	Modulation Pulse Width & Pulse Rate
P5	0~50 mA	150-260 uS	100 Hz	Modulation Pulse Width
P6	0~50 mA	260 uS	7-80 Hz	Modulation Pulse Rate
P7	0~50 mA	150-260uS	120 Hz	Modulation Pulse Width
P8	Repeat Cycle P1-P7			

**SYSTEMS INCLUDE:**

Stimulator, lumbar belt (fits 40", up to 52" with expander), 1 set of spring coil built-in leadwire cables, 1 pkg of 4 gel pads, 2 AAA or 9V battery, instruction manual & carrying case.

**MEDI-STIM** INC.

217 Industrial Court Wabasha, MN 55981 USA

1-800-363-STIM | Fax: 1-651-565-2410 | www.medi-stim.com